

uDirect™ Bluetooth® Pairing Guides



Pairing and connecting uDirect™ with your LG® cellular phone

Charge and turn ON both your cellular phone and uDirect™. Place them next to one another.

Turning on LG® Bluetooth®

1. Find the connectivity settings of your phone. With most LG models one of the following navigation options applies:

Menu > Settings & Tools > Bluetooth, or

Menu > Settings > Connectivity > Bluetooth > Paired devices, or

Menu > Tools > Bluetooth

Setting uDirect™ into pairing mode

To initiate the pairing process on uDirect, set uDirect into pairing mode by pressing and holding the Bluetooth button for about two seconds until the Bluetooth LED indicator is blue and flashing quickly. uDirect is now in pairing mode for two minutes or until the pairing is complete.

Pairing uDirect™ with your cellular phone


1. While in the Bluetooth menu, select one of the following options. Depending on the model, pairing selection may be:

Add new device, Search, or

Add New > Search

The phone will start searching for Bluetooth-enabled devices within range.

Note: If you have not activated the Bluetooth functionality, you may be asked if you wish to turn on Bluetooth. This will start the search process for Bluetooth devices. Your cellular phone should present a list of Bluetooth devices discovered.

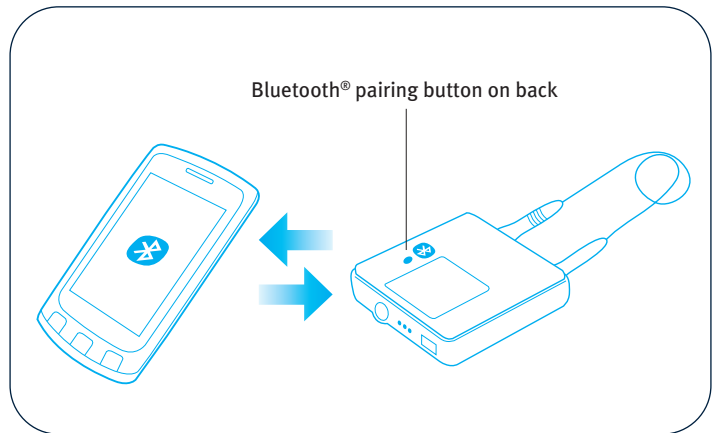

2. Under the list of available devices, select uDirect and enter 0000 (four zeros) when asked for the passcode or passkey. Allow devices to complete the pairing process. On the cellular phone, uDirect is shown as a paired device and on uDirect, the Bluetooth LED indicator flashes slowly, indicating that the devices are successfully paired.

Connecting uDirect™ with your cellular phone

In addition to being paired, which is done only once, your uDirect and cellular phone need to be connected. Typically, uDirect will be automatically connected after pairing.

1. To connect manually, in the Paired devices menu, select uDirect and

2. Select Options > Connect

On some phones you can set uDirect as authorized. This allows uDirect to connect automatically without confirmation.

3. If prompted to allow automatic reconnection, do so.

This way uDirect and your cellular phone will reconnect whenever within range.

If you have difficulties navigating to the Bluetooth settings, consult your cellular phone User Guide or ask your cellular phone dealer.