

Next™ 8 Custom

AutoPro2™


8 Channels, 8 Bands, Adaptive Directionality

HEARING INSTRUMENT FEATURES

- AutoPro2™ intelligently analyzes the input signal and quickly adapts to 1 of 2 distinct destinations. Within each destination, the adaptive feature can be customized for optimal listening and comfort
- Highly advanced feedback management that delivers more useable gain, allowing clients to enjoy the natural comforts and advantages of an open fit
- AntiShock™ instantaneously reduces the level of impulse noises such as a door slam, while maintaining the quality and intelligibility of speech
- Speech enhancement LD emphasizes speech signals based on the input level
- 8 channels provide high resolution signal processing
- Adaptive directional microphone system tracks and suppresses moving noise sources, while focusing on sounds from the front
- Noise Reduction, Wind Noise Manager
- Data logging accurately records data on time spent in each program and listening destination. Volume control changes are also logged in manual and automatic programs.
- MyMusic™ enhances the music listening experience by bringing out the rich, full tones of music
- OnBoard™ control is easily configured as a volume control or program button
- Easy-t provides automatic switching to a dedicated telephone program
- Ideal volume indicator provides a beep notification when preferred gain is reached on the volume control
- Up to 3 additional manual programs provide customization for individual needs and preferences
- Low battery warning
- Start up delay
- On/Off by opening or closing the battery door
- Can be programmed using NOAH-compatible U:fit™ and Standalone U:fit fitting software v1.4 or higher
- Choice of processing strategies, WDRC or Linear Limiting

OPTIONS AND ACCESSORIES


- Remote control with volume control, program change button, and more
- Telecoil (T) or Microphone/Telecoil (MT) option can be set as one of the three manual programs

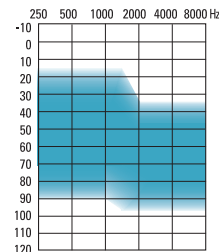

Fitting Guide




122/60
Full Shell Power

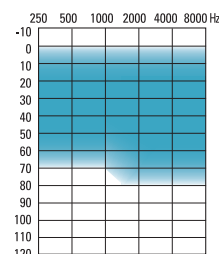

Fitting Guide




115/50
Full Shell

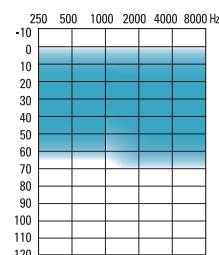

Fitting Guide



113/48
Half Shell / Canal

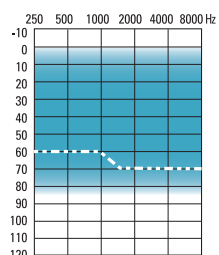

Fitting Guide



112/40
Mini Canal / CIC
116/55
CIC Power

--- Mini Canal / CIC
■ Power CIC

Next 8 Custom is suitable for fitting mild to severe hearing losses and can fit audiogram configurations ranging from reverse to precipitously sloping.

