

Stride™ V-PR Powered by Vivante




Stride V-PR

Stride V9-PR, Stride V7-PR, Stride V5-PR, Stride V3-PR, Stride V1-PR
Lithium-ion behind-the-ear (BTE) direct connectivity hearing instrument series

Integra OS		Level 9 Premium	Level 7 Advanced	Level 5 Standard	Level 3 Essential	Level 1 Lower Essential
Environmental classification	Total listening environments	8	6	4	2	AutoMic
	Conversation in loud noise	•				
	Conversation in a crowd	•				
	Conversation in a small group	•	•			
	Music	•	•			
	Noise	•	•	•		
	Quiet	•	•	•		
	Conversation in quiet	•	•	•	•	
	Conversation in noise	•	•	•	•	
	Total streaming environments	2	2	2	2	2
	Media speech	•	•	•	•	•
	Media music	•	•	•	•	•
	Sound optimization and performance	HyperFocus	•	•		
Conversation in loud noise (manual program)		•	•			
360 conversation in car (manual program)		•				
AutoFocus 360		•				
AutoFocus			•			
Dynamic noise reduction		•	•			
Soft speech lift		•	•	•		
Pinna Effect 2		•	•	•		
Pinna Effect		•	•	•	•	•
Speech enhancement		•	•	•	•	•
Noise reduction	•	•	•	•	•	
Hardware features	Tap control	•	•	•	•	
	Made for all direct connectivity	•	•	•	•	•
	Rechargeable	•	•	•	•	•
Fine-tuning channels	20	20	16	12	12	

Fitting guides

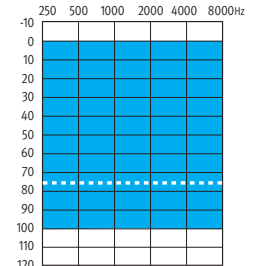

--- Slim tube

Available in all technology levels

Sound stabilization

- › AntiShock 2
- › Wind control
- › Feedback manager

Ease and convenience

- › Remote Plus app
- › TV Connector
- › Remote Control
- › PartnerMic
- › RogerDirect
- › DuoLink
- › Binaural Phone*

Microphone options

- › Adaptive directional
- › Fixed directional
- › Fixed wide directional
- › Omni directional

Fitting

- › First fit workflow
- › Automatic Adaptation Manager
- › Frequency compression 2
- › Tinnitus masker
- › Manual programs
- › IntelliVent

Experience innovations

Performance tracking

- › FLEX:TRIAL
- › Log It All
- › Ratings

Continuous care

- › Coach
- › Remote adjust
- › FLEX:UPGRADE

Personalization


- › Clarity and comfort buttons*
- › Optional app programs*
- › Equalizer

Stride V-PR is rated IP 68

Not all technology levels are available in all markets.


*Not available for lower essential level products.

ANSI/ASA 3.22 2014 (R2020)/IEC 60118-0: 2022 2cc coupler technical data

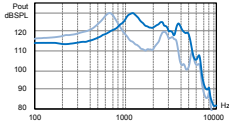


OSPL90

Maximum (dB SPL)	130	129
HFA - OSPL90 (dB SPL)	125	116

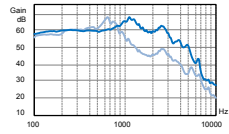

Full on gain (input 50 dB SPL)

Maximum (dB)	68	69
HFA - FOG (dB)	63	51

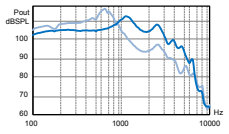

Basic Frequency Response

Reference test setting (RTS)

Frequency range (Hz)	<100 - 5700	<100 - 7100
Reference test gain (dB)	48	39
Typical battery life (h) *	18	18
Equivalent input noise at RTS (dB SPL)	19	19
Total harmonic distortion at 500 Hz/800 Hz/1600 Hz/3200 Hz (%)	4.0/3.0/2.0/1.0	1.0/1.0/2.0/1.0

Legend

- Earhook
- - - Slim Tube

General Test Information

Lithium-Ion rechargeable battery; Source: voltage 3.8 V
 * Typical operating time of the rechargeable battery is based upon a combination of bluetooth streaming and regular hearing instrument usage.
 The measurements obtained with a closed configuration using an HA-2 coupler (ANSI-3.7-1995).
 The hearing instrument set to Unitron TrueFit test settings. LLE is applied at an approximate level of 35 dB SPL.
 Domes should never be fit on patients with perforated eardrums, exposed middle ear cavities, or surgically altered ear canals. In the case of such a condition, we recommend use of a customized earmold.
 Monaural Latency in a fitted user mode is 6.5 mS according to ANSI 2051: 2017.
 We reserve the right to change specification data without notice as improvements are introduced.

WARNING: Changes or modifications to the hearing aid that are not explicitly approved by the manufacturer are not permitted. Such changes may damage the ear or the hearing aid.



Sonova AG
Laubisrütistrasse 28
CH-8712 Stäfa, Switzerland